

RT² Profiler PCR Array (Rotor-Gene® Format)

Mouse mTOR Signaling

Cat. no. 330231 PAMM-098ZR

For pathway expression analysis

Format	For use with the following real-time cyclers
RT ² Profiler PCR Array, Format R	Rotor-Gene Q, other Rotor-Gene cyclers

Description

The Mouse mTOR Signaling RT² Profiler PCR Array profiles the expression of 84 key genes involved in the mammalian target of rapamycin (mTOR) signaling pathway. mTOR, a serine/threonine protein kinase, integrates responses from a wide variety of signals (nutrients, hormones, growth factors and cellular stresses) to regulate cell growth, metabolism and survival. The first generation of mTOR inhibitors (e.g. rapamycin) failed to inhibit all mTOR functions, because the kinase forms two distinct protein complexes, mTORC1 and mTORC2. The rapamycin-sensitive mTORC1 complex regulates multiple biosynthetic cellular processes (protein synthesis, cell cycle progression, cell growth and proliferation). Until recently, the lack of mTORC2-specific inhibitors complicated elucidation of this protein complex's molecular functions. One definitive mTORC2 response is AKT activation, important for cell proliferation, migration and survival (apoptosis and autophagy inhibition). This array includes members of the mTORC1 and mTORC2 complexes as well as upstream regulators of many mTOR responses, and downstream genes from the many cellular processes regulated by mTOR complex activation. Using real-time PCR, you can easily and reliably analyze the expression of a focused panel of genes involved in mTOR signaling with this array.

For further details, consult the *RT² Profiler PCR Array Handbook*.


Shipping and storage

RT² Profiler PCR Arrays in the Rotor-Gene format are shipped at ambient temperature, on dry ice, or blue ice packs depending on destination and accompanying products.

For long term storage, keep plates at –20°C.

Note: Ensure that you have the correct RT² Profiler PCR Array format for your real-time cycler (see table above).

Note: Open the package and store the products appropriately immediately on receipt.

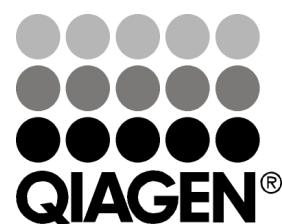

Sample & Assay Technologies

Array layout

The 96 real-time assays in the Rotor-Gene format are located in wells 1–96 of the Rotor-Disc™ (plate A1–A12=Rotor-Disc 1–12, plate B1–B12=Rotor-Disc 13–24, etc.). To maintain data analysis compatibility, wells 97–100 do not contain real-time assays but will contain master mix to account for weight balance.

Gene table: RT² Profiler PCR Array

Position	UniGene	GenBank	Symbol	Description
A01	Mm.6645	NM_009652	Akt1	Thymoma viral proto-oncogene 1
A02	Mm.148007	NM_026270	Akt1s1	AKT1 substrate 1 (proline-rich)
A03	Mm.177194	NM_007434	Akt2	Thymoma viral proto-oncogene 2
A04	Mm.235194	NM_011785	Akt3	Thymoma viral proto-oncogene 3
A05	Mm.26135	NM_133781	Cab39	Calcium binding protein 39
A06	Mm.179091	NM_026908	Cab39l	Calcium binding protein 39-like
A07	Mm.1022	NM_009861	Cdc42	Cell division cycle 42 homolog (S. cerevisiae)
A08	Mm.3996	NM_007700	Chuk	Conserved helix-loop-helix ubiquitous kinase
A09	Mm.21697	NM_029083	Ddit4	DNA-damage-inducible transcript 4
A10	Mm.250841	NM_030143	Ddit4l	DNA-damage-inducible transcript 4-like
A11	Mm.295397	NM_145470	Deptor	DEP domain containing MTOR-interacting protein
A12	Mm.290022	NM_145625	Eif4b	Eukaryotic translation initiation factor 4B
B01	Mm.3941	NM_007917	Eif4e	Eukaryotic translation initiation factor 4E
B02	Mm.6700	NM_007918	Eif4ebp1	Eukaryotic translation initiation factor 4E binding protein 1
B03	Mm.259516	NM_010124	Eif4ebp2	Eukaryotic translation initiation factor 4E binding protein 2
B04	Mm.278458	NM_008019	Fkbp1a	FK506 binding protein 1a
B05	Mm.141864	NM_010223	Fkbp8	FK506 binding protein 8
B06	Mm.394930	NM_019827	Gsk3b	Glycogen synthase kinase 3 beta
B07	Mm.3879	NM_010431	Hif1a	Hypoxia inducible factor 1, alpha subunit
B08	Mm.334313	NM_008284	Hras	Harvey rat sarcoma virus oncogene 1
B09	Mm.239865	NM_008300	Hspa4	Heat shock protein 4
B10	Mm.268521	NM_010512	Igf1	Insulin-like growth factor 1
B11	Mm.29254	NM_008343	Igfbp3	Insulin-like growth factor binding protein 3
B12	Mm.277886	NM_010546	Ikbbk	Inhibitor of kappaB kinase beta
C01	Mm.274846	NM_010562	Ilk	Integrin linked kinase
C02	Mm.4946	NM_008387	Ins2	Insulin II
C03	Mm.268003	NM_010568	Insr	Insulin receptor
C04	Mm.4952	NM_010570	Irs1	Insulin receptor substrate 1
C05	Mm.196581	NM_011949	Mapk1	Mitogen-activated protein kinase 1
C06	Mm.8385	NM_011952	Mapk3	Mitogen-activated protein kinase 3
C07	Mm.270866	NM_177345	Mapkap1	Mitogen-activated protein kinase associated protein 1
C08	Mm.289516	NM_019988	Mlst8	MTOR associated protein, LST8 homolog (S. cerevisiae)
C09	Mm.21158	NM_020009	Mtor	Mechanistic target of rapamycin (serine/threonine kinase)
C10	Mm.234502	NM_008659	Myo1c	Myosin IC
C11	Mm.10504	NM_011062	Pdpk1	3-phosphoinositide dependent protein kinase 1
C12	Mm.194127	NM_181414	Pik3c3	Phosphoinositide-3-kinase, class 3
D01	Mm.260521	NM_008839	Pik3ca	Phosphatidylinositol 3-kinase, catalytic, alpha polypeptide
D02	Mm.213128	NM_029094	Pik3cb	Phosphatidylinositol 3-kinase, catalytic, beta polypeptide
D03	Mm.229108	NM_008840	Pik3cd	Phosphatidylinositol 3-kinase catalytic delta polypeptide
D04	Mm.101369	NM_020272	Pik3cg	Phosphoinositide-3-kinase, catalytic, gamma polypeptide
D05	Mm.212039	NM_008875	Pld1	Phospholipase D1
D06	Mm.260177	NM_008876	Pld2	Phospholipase D2
D07	Mm.260288	NM_019411	Ppp2ca	Protein phosphatase 2 (formerly 2A), catalytic subunit, alpha isoform
D08	Mm.26134	NM_028392	Ppp2r2b	Protein phosphatase 2 (formerly 2A), regulatory subunit B (PR 52), beta isoform
D09	Mm.275393	NM_138748	Ppp2r4	Protein phosphatase 2A, regulatory subunit B (PR 53)
D10	Mm.207004	NM_001013367	Prkaa1	Protein kinase, AMP-activated, alpha 1 catalytic subunit
D11	Mm.48638	NM_178143	Prkaa2	Protein kinase, AMP-activated, alpha 2 catalytic subunit
D12	Mm.458152	NM_031869	Prkab1	Protein kinase, AMP-activated, beta 1 non-catalytic subunit
E01	Mm.31175	NM_182997	Prkab2	Protein kinase, AMP-activated, beta 2 non-catalytic subunit
E02	Mm.6670	NM_016781	Prkag1	Protein kinase, AMP-activated, gamma 1 non-catalytic subunit
E03	Mm.33649	NM_145401	Prkag2	Protein kinase, AMP-activated, gamma 2 non-catalytic subunit
E04	Mm.166501	NM_153744	Prkag3	Protein kinase, AMP-activated, gamma 3 non-catalytic subunit
E05	Mm.222178	NM_011101	Prkca	Protein kinase C, alpha
E06	Mm.207496	NM_008855	Prkcb	Protein kinase C, beta
E07	Mm.7980	NM_011102	Prkcc	Protein kinase C, gamma
E08	Mm.24614	NM_011104	Prkce	Protein kinase C, epsilon
E09	Mm.245395	NM_008960	Pten	Phosphatase and tensin homolog

Position	UniGene	GenBank	Symbol	Description
E10	Mm.319175	NM_053075	Rheb	Ras homolog enriched in brain
E11	Mm.757	NM_016802	Rhoa	Ras homolog gene family, member A
E12	Mm.275811	NM_030168	Rictor	RPTOR independent companion of MTOR, complex 2
F01	Mm.379007	NM_009096	Rps6	Ribosomal protein S6
F02	Mm.301827	NM_009097	Rps6ka1	Ribosomal protein S6 kinase polypeptide 1
F03	Mm.268383	NM_011299	Rps6ka2	Ribosomal protein S6 kinase, polypeptide 2
F04	Mm.220417	NM_153587	Rps6ka5	Ribosomal protein S6 kinase, polypeptide 5
F05	Mm.394280	NM_028259	Rps6kb1	Ribosomal protein S6 kinase, polypeptide 1
F06	Mm.271937	NM_021485	Rps6kb2	Ribosomal protein S6 kinase, polypeptide 2
F07	Mm.209933	NM_028898	Rptor	Regulatory associated protein of MTOR, complex 1
F08	Mm.31178	NM_178376	Rraga	Ras-related GTP binding A
F09	Mm.190922	NM_001004154	Rragb	Ras-related GTP binding B
F10	Mm.220922	NM_017475	Rragc	Ras-related GTP binding C
F11	Mm.300814	NM_027491	Rragd	Ras-related GTP binding D
F12	Mm.28405	NM_011361	Sgk1	Serum/glucocorticoid regulated kinase 1
G01	Mm.44231	NM_011492	Stk11	Serine/threonine kinase 11
G02	Mm.286006	NM_172656	Stradb	STE20-related kinase adaptor beta
G03	Mm.288702	NM_027880	Telo2	TEL2, telomere maintenance 2, homolog (S. cerevisiae)
G04	Mm.222	NM_011640	Trp53	Transformation related protein 53
G05	Mm.224354	NM_022887	Tsc1	Tuberous sclerosis 1
G06	Mm.30435	NM_011647	Tsc2	Tuberous sclerosis 2
G07	Mm.271898	NM_009469	Ulk1	Unc-51 like kinase 1 (C. elegans)
G08	Mm.162025	NM_013881	Ulk2	Unc-51 like kinase 2 (C. elegans)
G09	Mm.282184	NM_009505	Vegfa	Vascular endothelial growth factor A
G10	Mm.15607	NM_011697	Vegfb	Vascular endothelial growth factor B
G11	Mm.1402	NM_009506	Vegfc	Vascular endothelial growth factor C
G12	Mm.289630	NM_011739	Ywhaq	Tyrosine 3-monooxygenase/tryptophan 5-monooxygenase activation protein, theta polypeptide
H01	Mm.328431	NM_007393	Actb	Actin, beta
H02	Mm.163	NM_009735	B2m	Beta-2 microglobulin
H03	Mm.343110	NM_008084	Gapdh	Glyceraldehyde-3-phosphate dehydrogenase
H04	Mm.3317	NM_010368	Gusb	Glucuronidase, beta
H05	Mm.2180	NM_008302	Hsp90ab1	Heat shock protein 90 alpha (cytosolic), class B member 1
H06	N/A	SA_00106	MGDC	Mouse Genomic DNA Contamination
H07	N/A	SA_00104	RTC	Reverse Transcription Control
H08	N/A	SA_00104	RTC	Reverse Transcription Control
H09	N/A	SA_00104	RTC	Reverse Transcription Control
H10	N/A	SA_00103	PPC	Positive PCR Control
H11	N/A	SA_00103	PPC	Positive PCR Control
H12	N/A	SA_00103	PPC	Positive PCR Control

Related products

For optimal performance, RT² Profiler PCR Arrays should be used together with the RT² First Strand Kit for cDNA synthesis and RT² SYBR[®] Green qPCR Mastermixes for PCR.

Product	Contents	Cat. no.
RT ² First Strand Kit (12)	Enzymes and reagents for cDNA synthesis	330401
RT ² SYBR Green ROX [™] FAST Mastermix (2)*	For 2 x 96 assays in 96-well plates; suitable for use with the Rotor-Gene Q and other Rotor-Gene cyclers	330620

* Larger kit sizes available; please inquire.

RT² Profiler PCR Array products are intended for molecular biology applications. These products are not intended for the diagnosis, prevention, or treatment of a disease.

For up-to-date licensing information and product-specific disclaimers, see the respective QIAGEN kit handbook or user manual. QIAGEN kit handbooks and user manuals are available at www.qiagen.com or can be requested from QIAGEN Technical Services or your local distributor.

Trademarks: QIAGEN®, Rotor-Gene®, Rotor-Disc™ (QIAGEN Group); ROX™ (Applied Biosystems Corporation or its subsidiaries); SYBR® (Molecular Probes, Inc.).

1067688 03/2011 © 2011 QIAGEN, all rights reserved.

www.qiagen.com

Australia ■ 1-800-243-800

Austria ■ 0800/281010

Belgium ■ 0800-79612

Brazil ■ 0800-557779

Canada ■ 800-572-9613

China ■ 8621-3865-3865

Denmark ■ 80-885945

Finland ■ 0800-914416

France ■ 01-60-920-930

Germany ■ 02103-29-12000

Hong Kong ■ 800 933 965

Ireland ■ 1800 555 049

Italy ■ 800-787980

Japan ■ 03-6890-7300

Korea (South) ■ 080-000-7145

Luxembourg ■ 8002 2076

Mexico ■ 01-800-7742-436

The Netherlands ■ 0800 0229592

Norway ■ 800-18859

Singapore ■ 1800-742-4368

Spain ■ 91-630-7050

Sweden ■ 020-790282

Switzerland ■ 055-254-22-11

UK ■ 01293-422-911

USA ■ 800-426-8157



Sample & Assay Technologies