


Validation Data Sheet

QuantiNova LNA PCR Assay	MK167
Assay Catalog No.	SBH1220213
dPCR Instrument	QIAcuity digital PCR instrument
dPCR Master Mix	QIAcuity EG PCR Master Mix
Experimental Sample(s)	gBlock

Gene Expression Analysis: Concentration values [copies/ μ l] for different gBlock sample amounts (0 copies, 9000copies, 13500 copies and 18000 copies). Technical triplicates are shown.

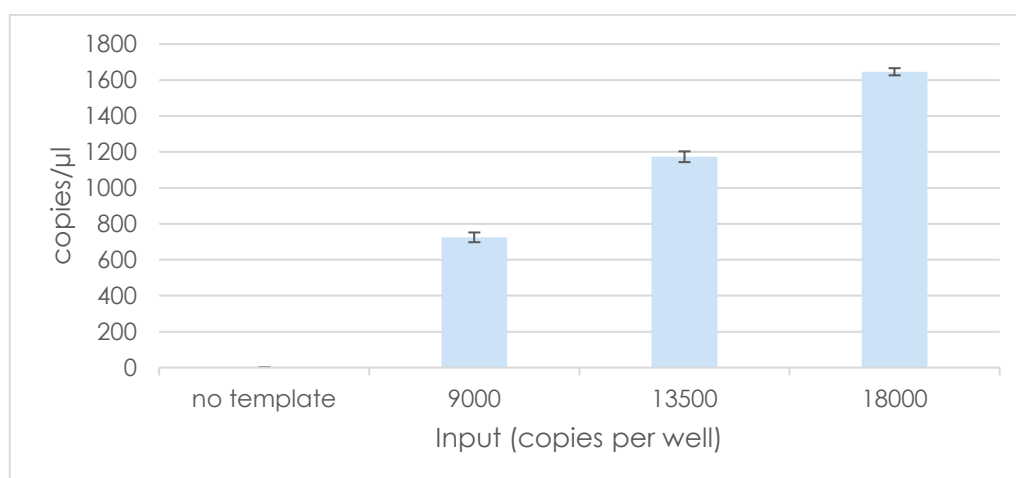

1D Scatter Plot: Single-well data (representative well of technical triplicates) for the different gBlock input amounts (FAM channel for EvaGreen detection; gBlock amounts in the 4 wells: 0 copies/ μ l, 750 copies/ μ l, 1125 copies/ μ l and 1500 copies/ μ l).

